

The ARS **PS-L Probe Station** is designed for the ultimate in flexibility for non-destructive device testing. The ability to upgrade and modify this system for further device testing has been in the forefront of this design.

The core of the system is the same for both the Liquid and Closed Cycle Cryogenic Probe Stations. These systems share the same table top, vacuum chamber, and probe arms. This enables the Flow Cryostat to be removed and replaced with a closed cycle (cryogen free) configuration at a later point.

This model of Probe Station uses one of ARS's Liquid Helium/Nitrogen flow cryostats. This allows for sample stage temperatures of ~4K or ~77K depending upon the cryogen that is used.

This system is designed to provide a large, clean sample environment. The vacuum chamber is made out of welded stainless steel and the radiation shield is made out of nickel plated OFHC copper. The low emissivity of the nickel plating and the high conductivity of the copper allow for a colder radiation shield and larger net refrigeration at the sample space. The high quality vacuum components are critical as it allows for deeper vacuum levels and cleaner samples with better electrical contact.


ARS' integrated approach of manufacturing for both the crycooler and the probe station ensures consistent performance, and also facilitates diagnostics and service of the integrated system.

Applications

- Microwave Properties
- DC, RF Properties
- MFMS
- Nanscale Electronics
- Superconductivity
- Electrical and Optical Properties of Nano Circuits
- Quantum Dots and Wires
- Non Destructive Device Testing

Typical Configuration

- Liquid Flow Cryostat
- Transfer Line Standard 6' (8' or 10' optional)
- Flow Meter Panels
- 10" Diameter Stainless Steel Vacuum Chamber with 6 Micromanipulated Probe Ports and 2 Spare NW 80 Accessory Ports Mounted on an Anodized Aluminum Table Top with Extruded Aluminum Support System
- 8" Diameter Nickel Plated OFHC Copper Radiation Shield
- 2.25" Diameter Grounded Sample Holder made out of OFHC Copper
- DC, Microwave, or Fiber Optic Probes
- 4 Sensors and 2 Heaters for temperature control and monitoring
- Lakeshore LS336 (4) Channel Temperature Controller with interconnecting cable to the cryostat
- 7:1 Zoom Microscope with <2 micron resolution and coaxial or ring light. Includes a High Resolution 24" Wide Screen LCD Display and a light source for the microscope</p>

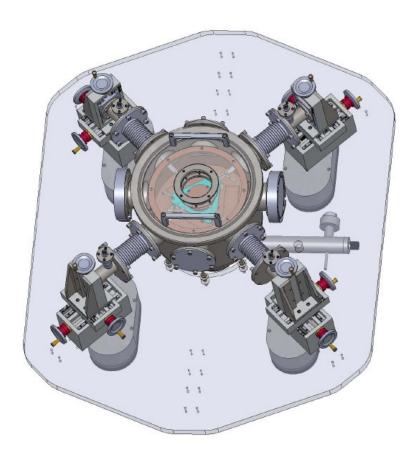



Cooling Technology

LT3	Open Cycle Cryostat
Liquid Cryogen Usage	Liquid Helium or Liquid Nitrogen

Temperature Range*

Liquid Helium	~4K - 400K **
Liquid Nitrogen	~77K - 400K **
Stability	0.1K (Can tune better with PID Control)

^{**500}K and 800K Options are also available

Radiation Shield

Material	Nickel Plated OFHC Copper
Diameter	8" (203 mm)
Ports	8 (standard) Empty ports have removable blanks installed
Lid	Removable Lid with Sapphire Cold Window
Mounted On	G-10 Pedestal
Optical Access	
Window Material	Sapphire
Window Diameter	2.5" (63 mm)
Window Clear View	2.25" (57 mm)

Sample Holders*


Grounded	2.25" Diameter (standard)
Electrically Isolated	2.25" Diameter (Optional)
Coaxial (Biased)	2.25" Diameter (Optional) Includes BNC Feedthrough with Coaxial Cable to sample stage
Triaxial (Biased with Guard)	2.25" Diameter (Optional) Includes Triaxial Feedthrough with Coaxial or Triaxial wire.
*Custom sample holders also available.	

Vacuum Chamber

Material	Welded, Stainless Steel
Diameter	10" (254 mm)
Probe Ports	6 (standard) 7 or 8 available upon request
Spare Accessory Ports	(2) NW-80 Flanges installed with blanks
Lid	Removable Lid with Quartz Window
Mounted On	29" x 36" Anodized Aluminum Plate
Optical Access	
Window Material	High Purity Quartz
Window Diameter	2.75" (69 mm)
Window Clear View	2.5" (63 mm)



The above picture shows the sample space of the probe station.

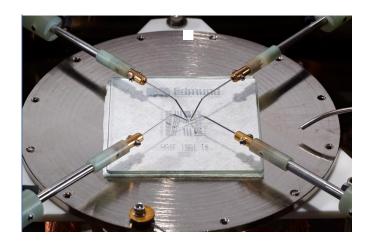

The above picture shows a custom sample holder and 4 DC probes.

^{*}Based on bare cold head with a closed radiation shield, and no additional sources of experimental or parasitic heat load



Translation Stages

Drive	Direct Drive
Bellows	Stainless Steel, Edge Welded
Motion	
X-Motion (Axial)	2" (50 mm) Standard
Y-Motion (Lateral)	1" (25 mm) Standard
	2" (50 mm) Optional
Z-Motion (Vertical)	0.5" (12.5 mm) Standard
Graduations	10 micron
Sensitivity	5 micron

Cool Down and Pump Down Time

Pump down Time* (<5 mTorr)	
Mechanical Pump (VPS-2)	~ 45 minutes
Turbo Pump (VPS-3)	~ 10 minutes
Cool down time to ~ 5 K**	~ 30-45 minutes

*The pump down time listed are approximate timelines using an ARS supplied vacuum pumping system. Actual vacuum levels of the system will be dependent of the vacuum pump itself. The system is capable of allowing for vacuum levels of 10^-6 Torr with an appropriate vacuum pumping system.

**The cool down time to 5K is listed for the cryocooler under typical conditions. If customizations are made, 5K sample stage temperatures may not be achievable.

Vibration Levels

the < 1 micron vibration levels.

Sample Stage Vibrations	< 1 micron
Vibration Dampening	Self Leveling Feet*
•	hen the feet are hooked up to a



Instrumentation for Temperature Control*

Temperature Sensors	Silicon Diodes
DT-670B-SD	Installed on the radiation shield
DT-670B-SD	Installed on the underside of the sample chuck for temperature control
DT-670B-SD	Installed on the cryocooler cold tip for diagnostics
DT-670-CU-4M	Calibrated Silicon Diode installed on top of the sample stage for accurate temperatures
Heaters	
50W Cartridge Heater	Installed on the under side of the sample stage for temperature control
(2) 50W Cartridge Heaters	Installed in parallel on the radiation shield to allow for 100W of heater power to assist in fast warm up of the system
Temperature Controller	LS-336 4 Channel Temperature Controller with 150W heater power

^{*}Interconnecting cables from the probe station to the temperature controller are supplied.



GSG Microwave Frequency Probe Arms*

*All GSG Microwave Frequency Probe Arms include the $\,$ Translation Stages. Probe arms are thermally anchored to the sample chuck and include +/- 5° Theta Planarization

40 GHz	Optional
Connector	K
Cable	Semirigid Coaxial
Frequency	0 to 40 GHz
Tip Material	Tungsten Beryllium Copper
Pitch	50 - 2450 micron (100 or 150 micron typical)
50 GHz	Optional
Connector	2.4
Cable	Semirigid Coaxial
Frequency	0 to 50 GHz
Tip Material	Tungsten Beryllium Copper
Pitch	50 - 1250 micron (100 or 150 micron typical)
67 GHz	Optional
Connector	1.85
Cable	Semirigid Coaxial
Frequency	0 to 67 GHz
Tip Material	Tungsten Beryllium Copper
Pitch	50 - 1250 micron 100 or 150 micron typical

Fiber Optic Probe Arms*

Range	UV/VIS or VIS/IR
Connector	Male SMA 905
Sample Termination	Bare Polished
Size (Typical)	100 micron - 400 micron
Mode	Single Mode or Multi Mode
Cable Material	Polyimide, Fused Silica
*All Fiber Probe Arms include the Translation Stages This does	

not include a light source or detector.

*All DC/Low Frequency Probe Arms include the Translation Stages. Probe arms are thermally anchored to the sample chuck.

Microminiature Coax Cable	Standard	
Connector	SMA or BNC	
Frequency	0 to 100 MHz	
Impedance	50 Ohm	
Includes Outer Ground Shield with Clip Connector		
Triaxial Cable	Optional	

metades outer Ground Smeta with etip connector		
Triaxial Cable	Optional	
Connector	Triaxial - 3 Lug	
Frequency	0 to 100 MHz	
Impedance	50 Ohm	
Kelvin Probes**	Optional	
Cable	Coaxial or Triaxial	
Connector	SMA, BNC, or Triaxial	
Frequency	0 to 100 MHz	
Probe Tips:		
Material	Tungsten (Standard) Gold Plated Tungsten (Optional) Beryllium Copper (Optional)	
Tip Radius	0.5 micron (standard) Other radius' also available	

^{**}Kelvin Probes use 2 cables and 2 connectors, but converge down to a single tip.


GigE Microscope with Zoom Lens*


•		
Zoom	7:1 (Standard)	16:1 (Optional)
Sensor	1/2" CMOS	1/2" CMOS
Field of View	4.2 mm - 0.61 mm	12.8 mm - 0.8 mm
Working Distance	89 mm	89 mm
N.A.:	0.024 - 0.08	0.0090 - 0.15
Light (Includes Light Source)	Ring or Coaxial	Ring of Coaxial
Resolution	3 microns	2 microns
Stand	Boom Stand with XYZ manipulation	Boom Stand with XYZ manipulation
Computer Interface	Ethernet Cable	Ethernet Cable
High Resolution Monitor	24"	24"

DC/Low Frequency Probe Arms*



PS-L Probe Station Outline Drawings

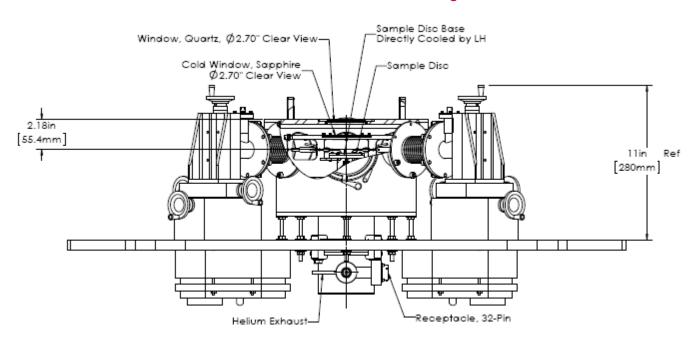

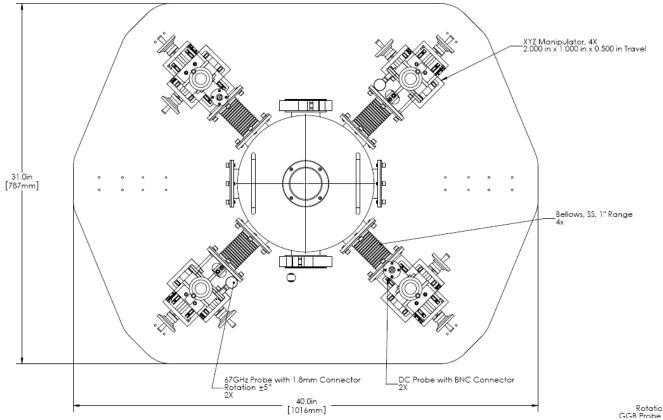

The above images shows a side profile of a probe station as well as a top down view.